
User Manual

Vitafit Smart Body Fat Scale

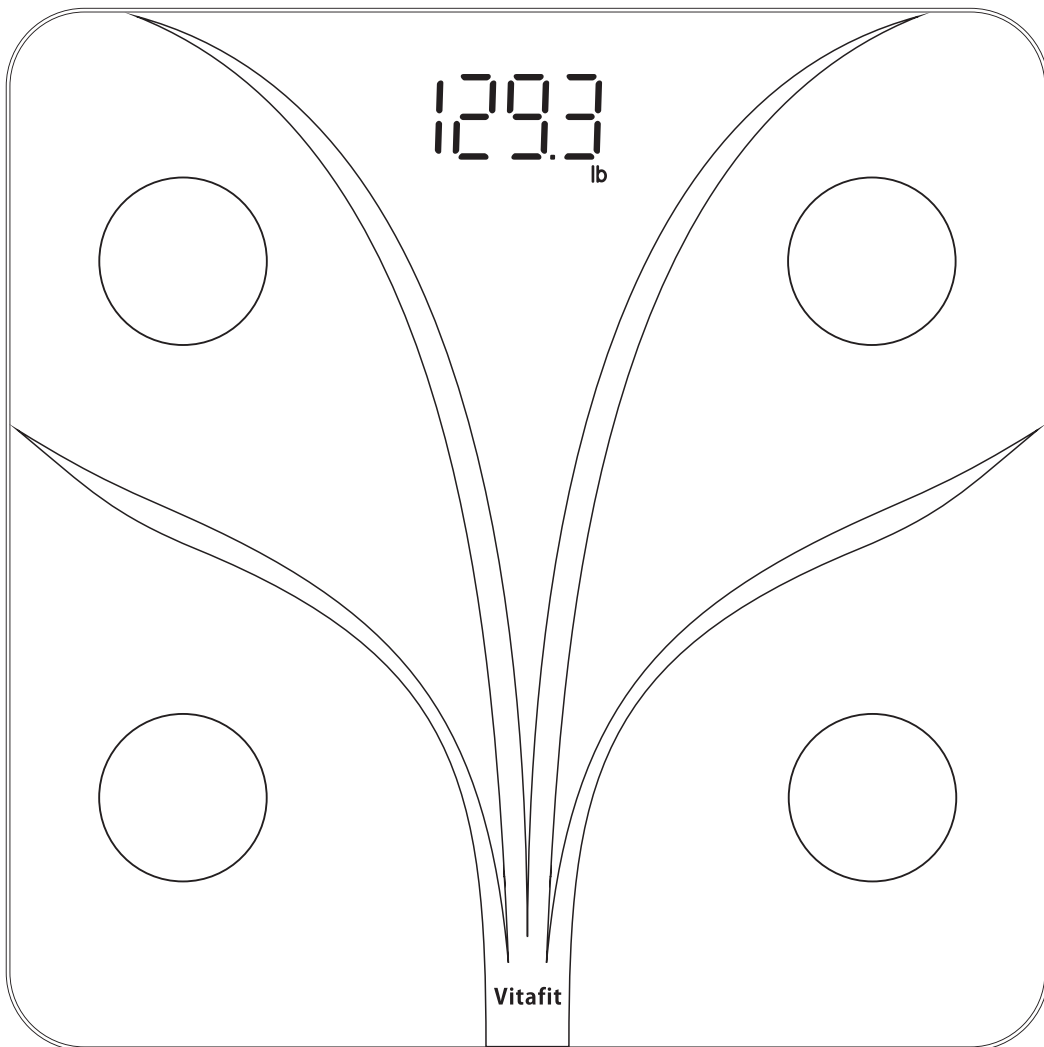Model:VT702

Please read this User Manual carefully before using the device

CONTENTS

01-05		English
06-10		Deutsch
11-15		Français
16-20		Español
21-25		Italiano

Warning

- ## Vitafit Smart Body Fat Scale

Safety Cautions

- ## Use and Care

- ## First Time Setup

[illegible]

Operation

Weighing without VitaFit App

- 1. Remove the isolator tab from the battery compartment.
- 2. For best accuracy, place the scale on a hard and flat floor. Carpeted or uneven floors can negatively effect the accuracy.
- 3. The scale operates with step-on technology but needs to be calibrated automatically under following occasions.
 - A. when first time use, moving the scale or battery exchange, the scale will display scale itself weight, and appear "CAL" after around 1second, and then display "0.0" or "0.00" after around 2 seconds , now it's ready to step on to weigh normally.

- B. If step on the scale before calibration, the scale will display an wrong weight and appear "CAL" after around 1second, and then display "0.0" or "0.00" after around 2seconds.

After calibration, please step again on to weigh and get correct weight.

- 4. Under other occasions, please step on to weigh directly.
- 5. Step off the scale, the weight readings remain for several seconds. Then the scale will turn off automatically.
- 6. When Weighing without App,the scale only display weight data and switch the unit by press the unit button on the backside of scale.

Weighing with VitaFit App


The VitaFit Bluetooth Scale can connect with your mobile phone or Ipad via Bluetooth . The Body Weight and body fat result data will be showed at the App named "VitaFit", which is easy to download from IOS App Store or Google Play.

Download App and Register ID

Scan the QR code or Download " VitaFit" via IOS App Store or Google Play.


Register Your own ID,Then add personal data and confirm OK.

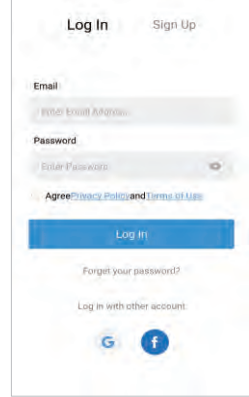Set Special User Mode


Under these modes, we adopt zero-current technology that only measure weight and BMI.

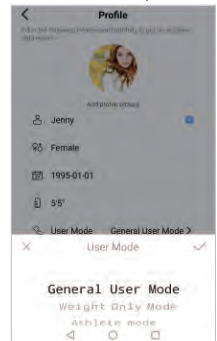- 1.If you are pregnant ,have a pacemaker or other internal device, please choose Weight Only Mode for safety measurement.

- 2.If the user is under 2years old, put in baby's age and height, VitaFit app will automatically recognize and enter baby mode that only measure weight.

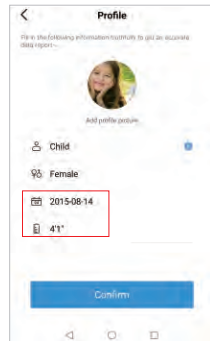- 3.If the user is among 2-18years old, put in child's age and height, VitaFit app will automatically recognize and enter child mode.

Note: When measuring with VitaFit App, please switch the unit via App,the scale will keep same unit with App. Let's make more clear,the unit button on the backside of scale will be nonfunctional when the scale connected to VitaFit App.

Pairing the Vitafit Scale with your smartphone


Turn on Bluetooth® on your smartphone.

NOTE: For Android users, Location and Bluetooth® must be turned on for Android 6.0 and later.

OR ELSE THE APP DON'T GET RESULTS ACCORDING TO ANDROID REGULATIONS.

1. Tap + on the measuring page to discover the scale.

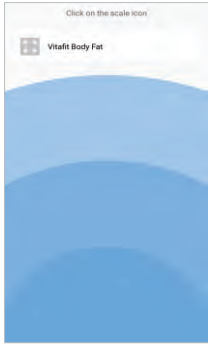2. Step onto the scale with bare feet, then add the device named "Vitafit Body Fat".

3. Keep standing on the scale for 5-10 seconds. Pairing is successful when your final weight and body fat result appear on the Vitafit app.


Viewing Reports

1. The body measurement result divided into two part "Body composition and Body Index", Click each bar to view each part of result.

2. Click each measurement result bar, you can view history data for the previous three times.

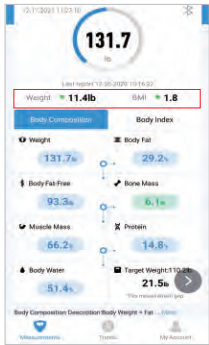
3. To view comparison data, click the bar under the weight section.

4. To view comparison data for a specific time, select starting date and ending date to generate the report, you can save or share your progress with friends via Facebook and Instagram by Click the "share" button on this page.

5. Click trend to view your progress history. You can select your progress for the week, month, or year.

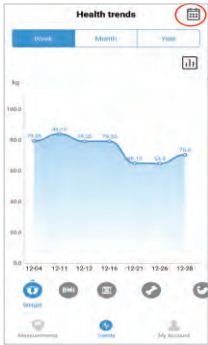
6. To view data for a specific day, click calendar logo the top right corner of the screen.


7. To view reference data from a specific time frame, click Data Select, then select date.

Deleting Reference Data

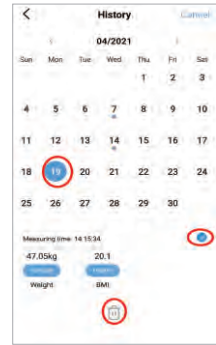
1. To delete a reference data, click Select and select date, then click the "X".





Adding/Deleting Users and Other Setting

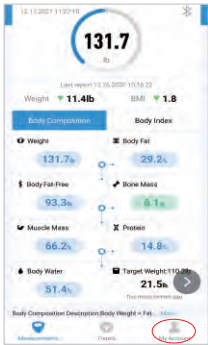1. On the Measurement page, tap My Account at the bottom of the screen.

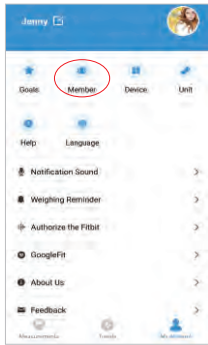2. To add new users, click Member and select "Add a user" and then enter the new user information.

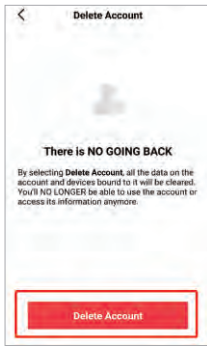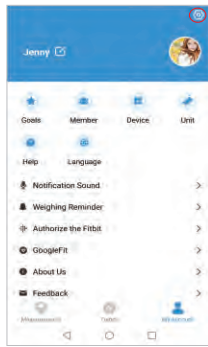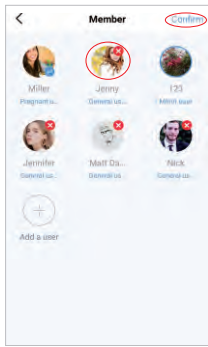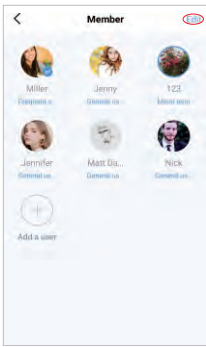3. To delete a user, first click Edit and delete, and then confirm.

4. To delete the account, first click SET logo and click" Delete Account" .




5.other setting: You can set goal weights, change weight units, select sounds, select languages, and set weighing reminder on the current page.

Connecting with Other Fitness Apps

Google Fit

1. On the page of My Account, Select Google Fit from the menu.
2. Turn on Google Fit, then sign into your Google Account. You will be automatically redirected to the Google Fit page. After successfully connecting, your measurement data will sync with Google Fit.




Note: In Android system, couldn't be allowed to Sync with Apple Health.

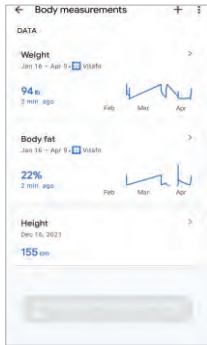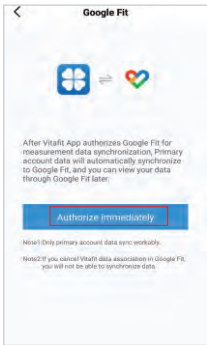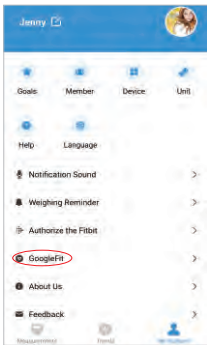Apple Health




1. On the page of My Account, Select Apple Health from the menu.
2. Click Authorize and the App will authorize immediately.
3. Open the Health App on your iPhone, then tap Sources. Be sure to turn on your Bluetooth®.
4. Select Vitafit from the list.
5. Turn on all categories to allow the Vitafit app to work.

Note: In IOS system, couldn't be allowed to Sync with Google Fit.

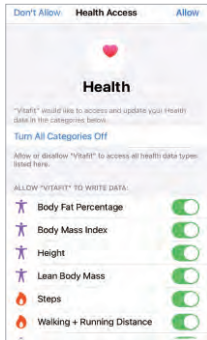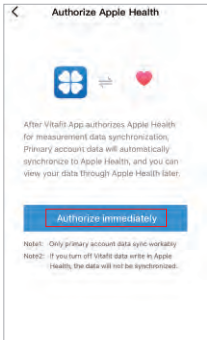
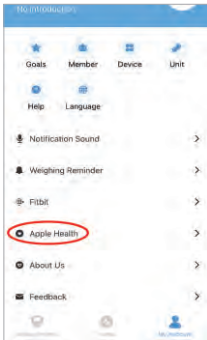
Fitbit

1. On the page of My Account, Select Fitbit from the menu.
2. Turn on Fitbit, then sign into your Fitbit account.

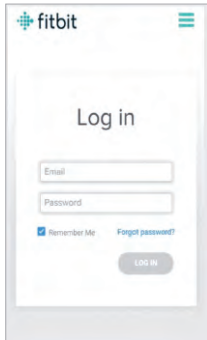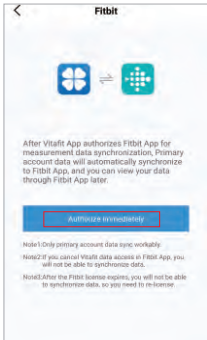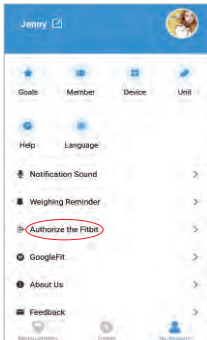3. Tap Allow to access the Fitbit app.
4. After successfully connecting, your data will sync to the Fitbit app. You will be able to view your progress both on the app as well as online.

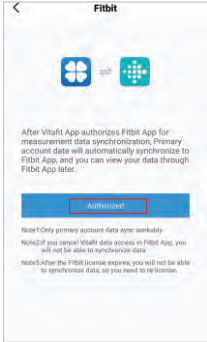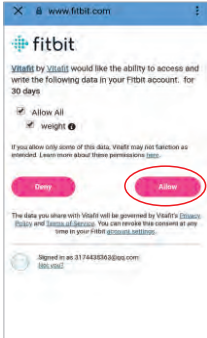Frequently Asked Questions

1. Why doesn't the scale turn on?
 - Check whether the batteries are properly installed.
2. Why does the data flash on the screen and then disappear?
 - The batteries are low and to be replaced.
3. Why doesn't the scale connect with the VitaFit app?
 - Make sure your phone's operating system is running on iOS 9.0+ or Android 4.3+.
 - Download and open the latest version of the VitaFit app.
 - Turn on Bluetooth® on your phone (for Android users, make sure your Location is turned on as well).
4. Why is my body fat not measured when weighing in?
 - Step on the scale with dry and bare feet, Place bare feet on the scale's conductive areas.
 - Make sure Bluetooth® is turned on.
 - Enter your personal information when creating your VitaFit account.
 - Make sure to place your feet on the scale's conductive areas.

Package Contents

Serial Number	Item	Quantity
1	Smart Body Fat Scale	1
2	1.5V AAA Battery	3
3	User Manual	1

Specifications

Capacity: Min 5 kg/11lb/0.11lb
 Max 182kg/400lb/28st

Weight unit: kg/lb/st

Graduation: 0.05kg/0.1lb/0.1lb

Dimensions: 280x280x24mm or 300x300x24mm(Depended on the model)

Display: LCD or LED (Depended on the model)

Platform: Tempered glass

Batteries: 3xAAA

Warranty Information

VitaFit guarantees all products to be of the highest quality in material and craftsmanship, and provides all products with 3 years warranty effective from the date of purchase.

For more information please visit www.51vitafit.com, www.vitafit.com

Email: service@51vitafit.com

Notice:	
	This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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Changes or Modifications not expressly approved by the party responsible could void the user's authority to operate this device.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement.

This device complies with Industry Canada Licence-exempt RSSs. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Vitafit

Inspired by Life

Factory: SHENZHEN SHINE INDUSTRIAL CO., LTD.
Address: 2&3F, Bldg5, 1st Industrial Zone, Changzhen Community,
Yutang Street, Guangming District, Shenzhen City, Guangdong
Province, China, 518132
Tel: +86-755-27175093
E-mail: sales03@iascent.com.cn

Company: Apex CE Specialists GmbH
Address: Habichtweg 1 41468 Neuss Germany.
E-mail: info@apex-ce.com
Tel: +353212066339

Company: APEX CE SPECIALISTS LIMITED.
Address: 89 Princess Street, Manchester, M1 4HT, UK
E-mail: Info@apex-ce.com
Tel: +441616371080